


POLYMATIC DUAL HVE TIP

Provides One-Handed Aspiration Plus Tongue & Cheek Retraction.



Polymatic keeps the cheek and tongue out of the way reliably and provides a clear view of the procedure.

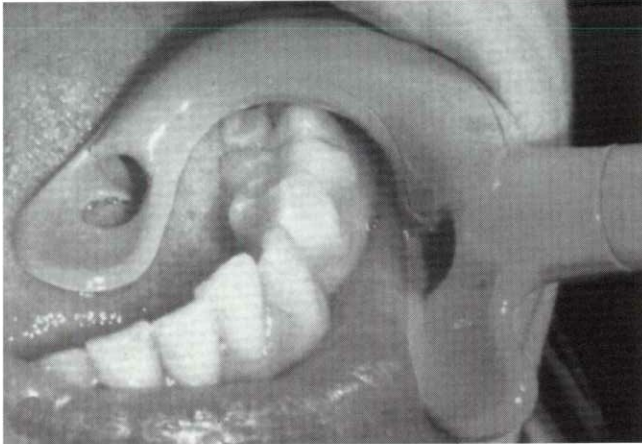

Polymatic's unique design is easy to handle in all types of treatment of the lower jaw, especially for the removal of calculus.

Aspirates the Buccal and Lingual at the Same Time.

- **Same Time Retraction of the Cheek and Tongue.**
- **No Aspiration of Soft Tissue.**
- **Can be Autoclaved.**
- **Comes in a Set of 2 for Left and Right Side Aspiration.**
- **Fits all HVE Handpieces.**

Polymatic Dual Function Aspirator Tips can be used for any type of procedure with high/low speed handpiece, air/water syringe or sonic scaling. It aspirates spray mist and patient secretions with ease from the bottom of the mouth. The cheek and tongue are kept out of the way with a very simple one-handed holding technique of the dentist, assistant or hygienist. A clear view of the working area from the posterior/wisdom to the anterior/incisal teeth occurs. Suck down of the soft tissue around the tip opening is not possible.

Polymatic is made of a durable plastic adding a warming thermal buffer to the patient; unlike metal instruments which are cold. Polymatic is completely autoclavable up to 340°F (170°C), including the very demanding dry heat chambers.